




SWIFT

OPTICAL INSTRUMENTS

Exceptional Optics ■ Durable Construction ■ Innovative Designs

Wired Digital Microscope or Camera Running Simultaneously on Wired Network & Wireless Router

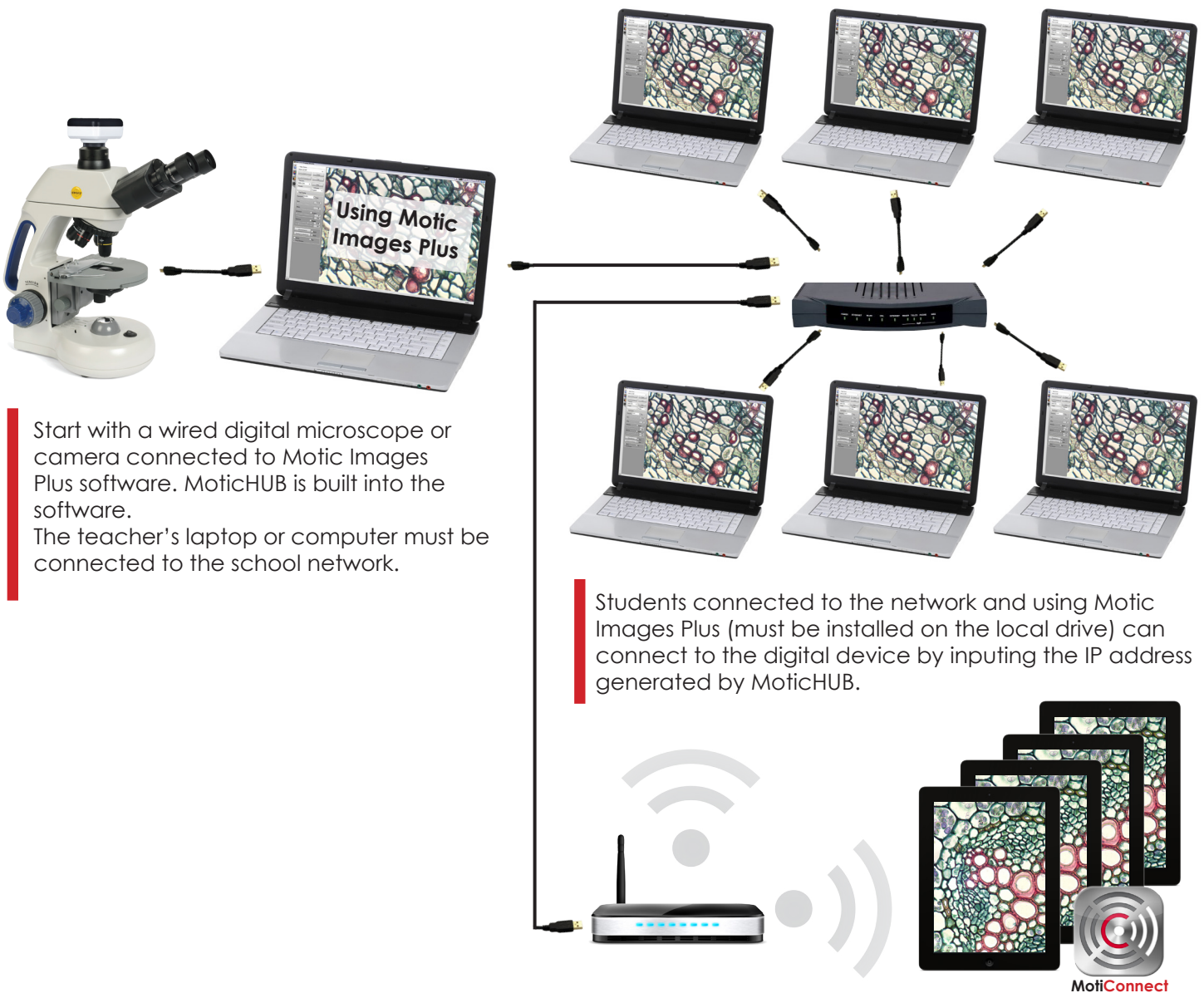

Start with a wired digital microscope or camera connected to MotiC Images Plus software. MotiHUB is built into the software. The teacher's laptop or computer must be connected to the school network.

Students connected to the network and using MotiC Images Plus (must be installed on the local drive) can connect to the digital device by inputting the IP address generated by MotiHUB.

Students using tablets or smart phones can connect to the digital microscope or camera using the free MotiConnect app and the IP address generated by MotiHUB.